


UNIVERSAL RADIOLOGY SYSTEMS FOR OPTIMAL COST EFFICIENCY








Manufactured in an ISO/EN factory, our product range offers reliability without compromise.

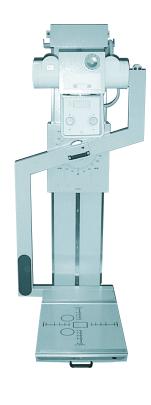

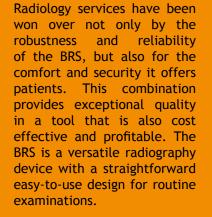

The FERIA consists of an U arm with multiple projections, perfectly counterbalanced, installed on a floor mounted column. Total patient coverage can be achieved thanks to the combination of the stand's movements as well as the wheels with brakes.

FEST

The FERIA is a universal radiography unit which allows for an important volume of exposures in all types of patient positions (standing, sitting or lying down).





With references on the 5 continents, the FERIA guarantees a cost effective system for a long term investment. The BRS responds in particular to the recommendations of the World Health Organisation (WHO).

FOR ALL TYPES OF EXAMS

Both the BRS and FERIA allow you to practice all types of exams in a standing, sitting or lying down position.

The following studies are available:

- → Lungs
- Bones
- General Radiography
- Skeleton
- Abdomen
- Extremities
- Emergency exams, etc ...



THE MAIN BENEFITS OF THE BRS / FERIA

- Oscillating anti-scattering grid (40l/cm ; R 10:1 according to detector)
- \bigcirc All cassettes from 13 x 18 cm to 36 x 43 cm are accepted
- Light diaphragm for optimal patient positionning
- →Motorized variable SID (99 to 183 cm) (FERIA)
- Tilting bucky (+/- 45°) (FERIA)
- \bigcirc Automatic Exposure Control (AEC) available as option
- These solutions meet the requirements of all types
 of radiology services, including those with heavy
 workflow.
- Thanks to its compactness, the FERIA's installation is very easy even in small rooms.
- Also available through our DR Retrofit Kit (EZ Pixium 3543) offering all the advantages of digital.











For more than 25 years, APELEM has been privileging the research and development of innovative X-Ray systems in order to supply our customers with equipment featuring cutting edge technology.

